

# AWS

## BLADE Scale User Manual



**BLADE-100**  
**BLADE-600**  
**BLADE-1KG**

## BLADE Scale


Thank you for purchasing the American Weigh BLADE digital pocket scale. Please read all operating instructions carefully before use. Scales are precision measuring devices and should always be handled with proper care. To ensure years of reliable service, keep these simple tips in mind:

- Do not exceed the scale's maximum capacity. Overloading your scale can permanently damage it!
- Avoid exposure to extreme heat or cold. Scales perform best at normal room temperature. If temperatures have changed dramatically, calibration and adjustment may be necessary.
- Allow your scale to warm up for 30-60 seconds before performing calibration and adjustment(if available).
- Store your scale in a clean, dry location. Dust, dirt, and moisture can accumulate on the weighing sensors and electronics causing inaccuracy or malfunction.
- Avoid using your scale in close proximity to strong radio wave emitting devices such as computers, cash registers, and cordless phones.
- Always weigh on a flat and level surface, free from vibrations and drafts. The corner of a room is usually the most stable.
- Gently apply all items to be weighed. Do not drop items onto the weighing platform.
- Avoid dropping your scale. The warranty does not cover damage due to rough treatment or overload.
- Check the batteries first if you are having any trouble with your scale. This simple step can remedy most scale issues.
- Do not disassemble your scale. This product contains no user serviceable parts.

## I. Part List



1. Scale
2. Cover/Expansion Tray
3. AAA Batteries (x2)

## II. Key Description

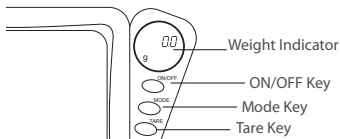

## III. Battery Installation

1. Before using the scale for the first time, check the battery compartment and remove any plastic insulation strips or packaging that may be in place to prevent battery drainage during shipping/storage.
2. If you are replacing the batteries, make sure the positive and negative contacts are properly aligned.

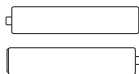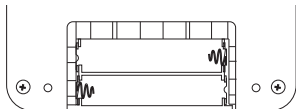

# BLADE Scale

## IV. Operation

1. Place the scale on a flat, level surface and press the **ON/OFF** key to turn the scale on.
2. Once the start-up process completes and the scale goes to zero you are ready to weigh.
3. Place the quantity to be weighed onto the weighing platform.
4. Wait a couple seconds for the display to stabilize before taking a weight reading.

## V. Changing the Unit of Measure

1. Press the **MODE** key to change the scales displayed unit of measure.
2. If the scale is turned off then back on, it will default to the last unit used.
  - **BLADE-100 Units:** g (gram) ▶ oz (ounce) ▶ ct (carat) ▶ gn (grain)
  - **BLADE-600, BLADE-1KG Units:** g (gram) ▶ oz (ounce) ▶ ozt (troy ounce) ▶ dwt (pennyweight)

## VI. Using the Tare Feature

You can use this scale's tare feature to subtract the weight of an empty container or vessel for net weight determination. **NOTE:** The weight of the container will subtract from the available capacity.

1. Place the empty container or vessel on the weighing platform.
2. Press the **TARE** key once. The display will then reset to zero. If the container is removed, it's weight will be displayed as a negative number.
3. You may now fill the container with the quantity to be weighed. The display will show the net weight of

the contents.

4. To return to normal weight display (gross weight), remove any items from the weighing platform and press the **TARE** key again to re-zero the scale.

## VII. Calibration

Calibration and adjustment are performed at the factory. However, the weighing range can shift slightly during shipping. Also, due to changes in local acceleration of gravity around the world, adjustment is recommended when moving the scale from one location to another. To check calibration:

1. Place the required calibration weight on the scale and note the displayed value.
2. If the displayed value is not within acceptable tolerance (see chart below), adjustment is required.

Model	Acceptable Tolerance ( $\pm$ )	Calibration Weight	Recommended Accuracy Class
BL-100	0.02g	100g	OIML M2
BL-600	0.2g	500g	OIML M2
BL-1KG	0.2g	500g	OIML M2

## VIII. Adjustment

1. Press and hold the **MODE** key until the display shows "CAL".
2. Press the **MODE** key again. The display will flash "CAL" followed by the calibration weight value (see chart).
3. Place the calibration weight on the center of the weighing platform. Make sure you are using the correct weight for your specific model.

## BLADE Scale

4. Once the weight is accepted, the display will show "PASS" then return to normal operation.

### IX. Troubleshooting

- **Problem:** Scale will not turn on
- **Solution:** Check the batteries
  
- **Problem:** Displayed weight fluctuates randomly
- **Solution:** Make sure your workspace is stable and free of any air currents or vibration
  
- **Problem:** Displayed weight is inaccurate
- **Solution:** Perform a calibration adjustment

### X. Error Codes

1. **OUTZ** - Zero range has shifted. Perform calibration adjustment. If the problem persists, the load cell may be permanently damaged due to overloading.
2. **O-Ld** - Maximum capacity exceeded. Remove the extra weight to avoid damaging the load cell.
3. **L** - Batteries are low. Replace the batteries.

## XI. Specifications

	BL-100
Max. Capacity	100g
Readability	0.01g
Min. Weight	0.03g
Power	2 x AAA
Dimensions	0.8 x 3.8 x 3.8"
Platform Dimensions	2.7 x 2.7"

	BL-600	BL-1KG
Max. Capacity	600g	1000g
Readability	0.1g	0.1g
Min. Weight	0.2g	0.2g
Power	2 x AAA	
Dimensions	0.8 x 3.8 x 3.8"	
Platform Dimensions	2.7 x 2.7"	

## XII. Capacity/Readability Charts

BL-100
100g x 0.01g
3.528oz x 0.001oz
500ct x 0.05ct
1543gn x 0.1gn

BL-600	BL-1KG
600g x 0.1g	1000g x 0.1g
21.16oz x 0.01oz	35.27oz x 0.01oz
19.29ozt x 0.01ozt	32.15ozt x 0.01ozt
385.8dwt x 0.1dwt	643.1dwt x 0.1dwt